

TABLE OF CONTENTS

Consumables	2
SIBO Testing	2
Equipment	8



EASYSAMPLER MOUTHPIECE/BAG ASM

EasySampler Mouthpiece/Bag ASM (must be purchased in Multiples of 50)

Used with the AlveoVac Extraction System.

SKU: QT.01665



BREATH COLLECTION VIALS

ALVEOVAC, COLLECTION VIALS, BOX OF 100.

Breath collection vials (must be purchased in Multiples of 10)

Used with the AlveoVac Extraction System.

SKU: QT.02602



ONE-WAY PLASTIC STOPCOCK

One-way stopcocks for breath testing syringes. (must be purchased in Multiples of 10)

Used with the BreathTracker.

SKU: QT.01727-V



EASYSAMPLER KIT

This test kit contains only a mouthpiece, 10 sample collection vials, 10 vial labels and general instructions for use. It does not contain any sugar or testing substrate.

Used with the BreathTracker.

SKU: QT.04207

FRUCTOSE, 25GM PACKET




Fructose Sugar 25gm Packet.
Used with the BreathTracker.
SKU: QT.02900-25-S

LACTEST ® - 25GM/PACKET



Unflavored Lactose Packets (25gm/ea)
Used with the BreathTracker.
SKU: QT.02425-S

SUGAR, SUCROSE, 50 GM PACKET

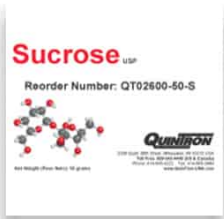


50 gram packet of unflavored Sucrose powder.
Used with the BreathTracker.
SKU: QT.02600-50-S

LACTULOSE, 10GM PACKET



Convenient, single-dose packets of cherry-flavored lactulose for Small Intestinal Bacterial Overgrowth (SIBO) and Intestinal Transit Time (ITT) breath tests. Used with the BreathTracker.
SKU: QT.02500-10-S

GAS COLLECTION BAG

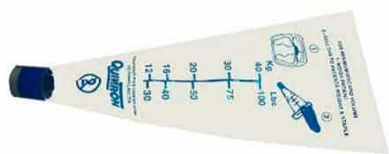

These bags are metallized to make it impermeable to gases, and laminated so it is heat-sealable. They have a one-way port to fit the Tee, and a Luer-fitting and stopcock for attaching a syringe. Bag, Collection, Single Patient, 750 ml (must be purchased in Multiples of 10). Used with the BreathTracker.

SKU: QT.00830-P

SAMPLE HOLDING BAG



These small 250mL bags are fabricated from the same gas impermeable material as the Collection Bag. Samples are transferred from the Collection Bag through the stopcocks (or by syringe), thus freeing the Collection Bag for additional sampling. 250 ml, Sample Holding Bag, Single Patient (must be Purchased in Multiples of 5) Used with the BreathTracker.

SKU: QT.00842-P


PLASTIC TEE-CONNECTOR



Plastic Tee-Connector, Single Patient (must be purchased in Multiples of 10). Used with the BreathTracker.

SKU: QT.00850-P

DISCARD BAG

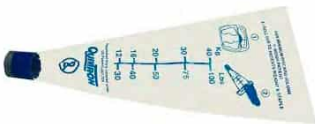

This disposable polyethylene bag is used to separate the “dead-space” air from alveolar air which goes into the Collection Bag.

Mouthpiece, 20mm, Blue, Moulded (must be purchased in Multiples of 10). Used with the BreathTracker.

SKU: QT.00843-P



DISPOSABLE MOUTHPIECES

Mouthpieces are inexpensive enough to be disposable, but they can be cleaned, disinfected and reused if desired.

Mouthpiece, 20mm, Blue, Moulded (must be purchased in Multiples of 10) Used with the BreathTracker.

SKU: QT.00991-P

QUINTRON BREATHTRACKER SC

The ability to analyse H₂, CH₄ and CO₂ provides physicians with necessary information to confidentially diagnose carbohydrate malabsorptions and/or small intestinal bacterial overgrowth (SIBO).

The BreathTracker SC analysers trace-gas concentrations of hydrogen (H₂), methane (CH₄) and carbon dioxide (CO₂) in a patients alveolar breath sample using solid-state sensors. The results are measured and displayed in parts per million (ppm) for H₂ and CH₄ and percent (%) for CO₂. The BreathTracker SC also utilise an added feature for sample contamination detection and correction based on CO₂ measured in patient samples.

SKU: QT.05019-M



INDEX

B

Breath Collection Vials 1

D

Discard Bag 6

Disposable Mouthpieces 7

E

EasySampler Kit 2

EasySampler Mouthpiece/Bag ASM 1

F

Fructose, 25gm Packet 3

G

Gas Collection Bag 5

L

LacTest ® - 25gm/packet 3

Lactulose, 10gm packet 4

O

One-Way Plastic Stopcock 2

P

Plastic Tee-Connector 6

Q

QuinTron BreathTracker SC 7

S

Sample Holding Bag 5

Sugar, Sucrose, 50 gm packet 4