
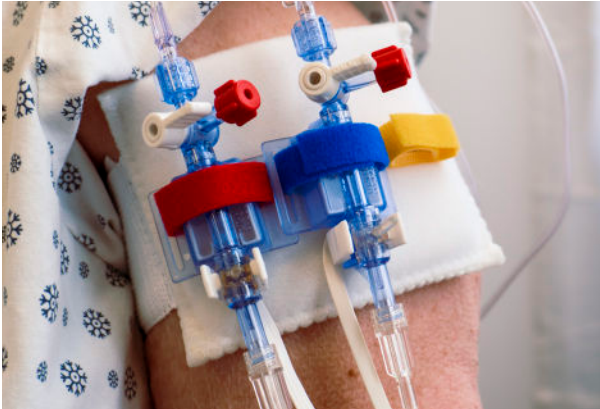


Transducer Holder

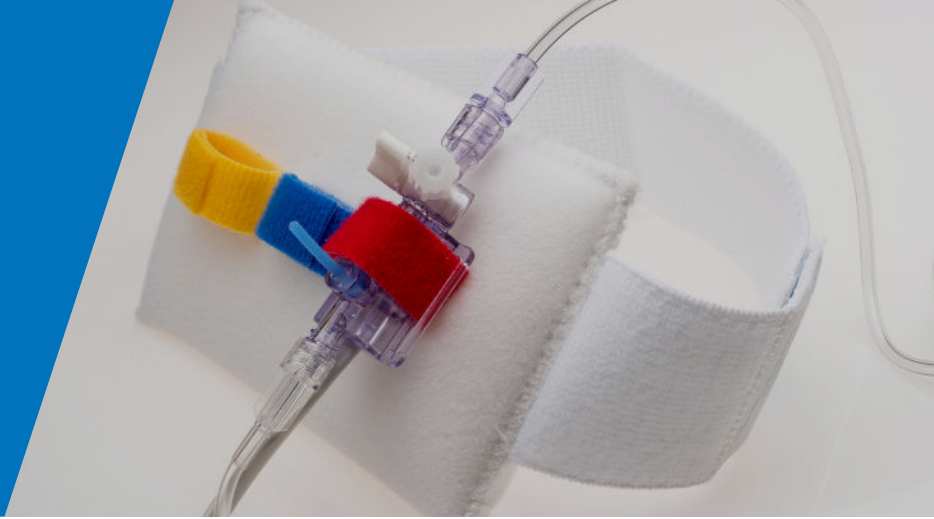

The Dale Transducer Holder reduces the likelihood of inaccurate measurements by maintaining proper positioning of each transducer at the phlebostatic axis. Made of a soft, comfortable material, the platform provides the base for the transducer tab holders and is secured to the patient's arm with an elasticized strap with adjustable hook and loop tabs. The Dale Transducer holder replaces the need to tape transducers on a patient's chest or mount transducers on an IV pole. Once in place, the Dale Transducer Holder eliminates the need for repeated leveling during patient repositioning, transport, or ambulation.

Key Features

- ◆ Stabilizes up to three transducers.
- ◆ Improves patient monitoring.
- ◆ Quick and easy to apply and remove.
- ◆ Saves nursing time.
- ◆ Can accommodate VAMP system (Venous Arterial blood Management Protection)
- ◆ Not made with natural rubber latex.

Product #	Description	Quantity
620	Fits up to three standard transducers One size fits most	10/box
622	Fits up to two large transducers One size fits most	10/box

Patents dalemed.com/patents

Transducer Holder

Features

Clinical Advantages

Secure Platform

- Secures and organizes up to three transducers, as well as the ability to accommodate VAMP system
- Platform is soft and padded
- Easy transducer positioning at phlebostatic axis enhances leveling accuracy
- Soft elasticized strap to secure to patient's upper arm without restriction
- Monitoring at phlebostatic axis is maintained during patient repositioning, transport, ambulation

Holding Tabs

- Application is easily identified, colors coordinate CVC, Arterial, and PA pressure lines
- Secures 1, 2, or 3 transducers
- Color coded to easily identify lines
- Easy to use hook and loop fasteners
- 620 - Secures up to 3 standard size transducers
- 622 - Secures up to 2 large size transducers
- 622 - Also compatible with a VAMP pressure line

Easy to Apply

- Reduces the time and trouble of releveling with patient or bed position changes
- Eliminates taping individual transducers to the patient's chest or arm or mounting on bedside pole
- Saves nursing time

Convenient System

- The holder is designed to secure and organize the transducers at phlebostatic axis
- Allows for a complete securement system for radial artery blood pressure monitoring when used in conjunction with the Dale Bendable ArmBoard

Not Made with Natural Rubber Latex

- Eliminates patient and caregiver latex reactions

Learn More

