

Freddy Therapy

Compact size, able to adapt to all situations



THERAPY CART

Frequent usage

Emergency room

Ambulatory

ICU

Dialysis

Hospitalization ward

Specialistic





DRAWERS IN STEEL SHEETS, SLIDING ON GLIDES WITH SOFT CLOSE LOCKING SYSTEM

Benefits



- Sturdy and with sheets in stainless steel painted with thermosetting paint.
- Accessories in stainless material.
- Single body structure and self-supporting.
- Thanks to its small size it is easy to handle even in narrow places.
- Despite its small size it is actually very capacious.
- Swivable wheels with sliding system for easy driving.
- The cart is upgradable with a wide range of accessories that can be placed on the DIN bars located on 3 sides.
- Push and pull extra tray.
- User friendly:
 - 3 different drawers configuration;
 - available up to 7 drawers;
 - the accessories can be changed or removed without any tools.
- Anti-tilting structure.
- Table top made in plastic material without sharp edges for easy cleaning and disinfection.
- Drawers in steel sheets, sliding on glides with soft close locking system
- Drawers centralized locking system.
- Glides can support up to 30 kg.

Standard configuration



- First drawer provided with guarantee seal.
- No. 10 self-destructing seals.
- No. 3 lateral DIN bars.
- No. 1 back DIN bar.
- Centralized locking system with key.
- Double key.
- No. 4 swivable and stainless wheels of which 2 with pedal brake. No marking, extremely silent, sliding on high performing roller bearing.
- No. 4 bumpers for each corner of the base.
- Push and pull extra tray.
- Table top with integrated push handle.



Do you want to customize the color?

Frame or click the QR code and find out how to customize the device by consulting our color card.



Recommended accessories

Kit composed by the following accessories:

- 1 Cod. ACCF1 - IV pole
- 2 Cod. ACCF7 - Bottle rack
- 3 Cod. ACCF8 - Case for gloves
- 4 Cod. ACCF9 - Needles rack
- 5 Cod. ACCF10 - Defibrillator or monitor tray



Recommended accessories



ACCF1



IV pole

ACCF2



Trash bin

ACCF3



Manual hand sanitizer

ACCF4



Automatic hand sanitizer

ACCF5



Hoxigen tank holder

ACCF6



Defibrillator or monitor trays



Recommended accessories



ACCF7



Bottle rack

ACCF8



Case for gloves

ACCF9



Needles rack

ACCF10



Defibrillator or monitor tray

ACCF12



Guarantee seal

ACCF13



Guarantee seal with serial number



Recommended accessories




ACCF14




Assemblable drawers cabinet

ACCF15

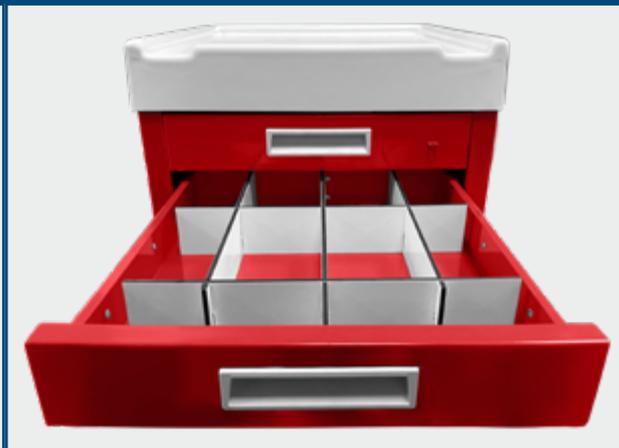

Drawer dividers: ACCF15/1 - 10 compartments; ACCF15/2 - 15 compartments; ACCF15/3 - 20 compartments; ACCF15/4 - 25 compartments

ACCF16

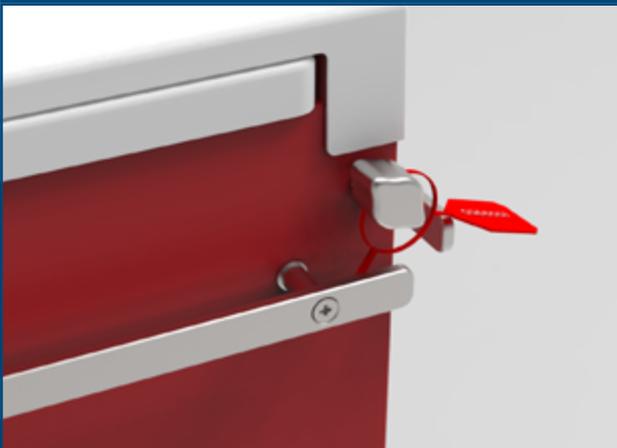

Centralized and sealable locking



Available drawer's configuration



3 drawers
dimensions 42x37,8x7,8h;

2 drawers
dimensions 42x37,8x15,8h;

1 drawer
dimension 42x37,8x23,7h;

GB0321.SP-1
Drawer's configuration 3+2+1



1 drawer
dimension 42x37,8x7,8h;

3 drawers
dimensions 42x37,8x15,8h;

1 drawer
dimension 42x37,8x23,7h;

GB0321.SP-2
Drawer's configuration 1+3+1



2 drawers
dimensions 42x37,8x7,8h;


1 drawer
dimension 42x37,8x15,8h;

2 drawers
dimensions 42x37,8x23,7h;

GB0321.SP-3
Drawer's configuration 2+1+2



Dimensions

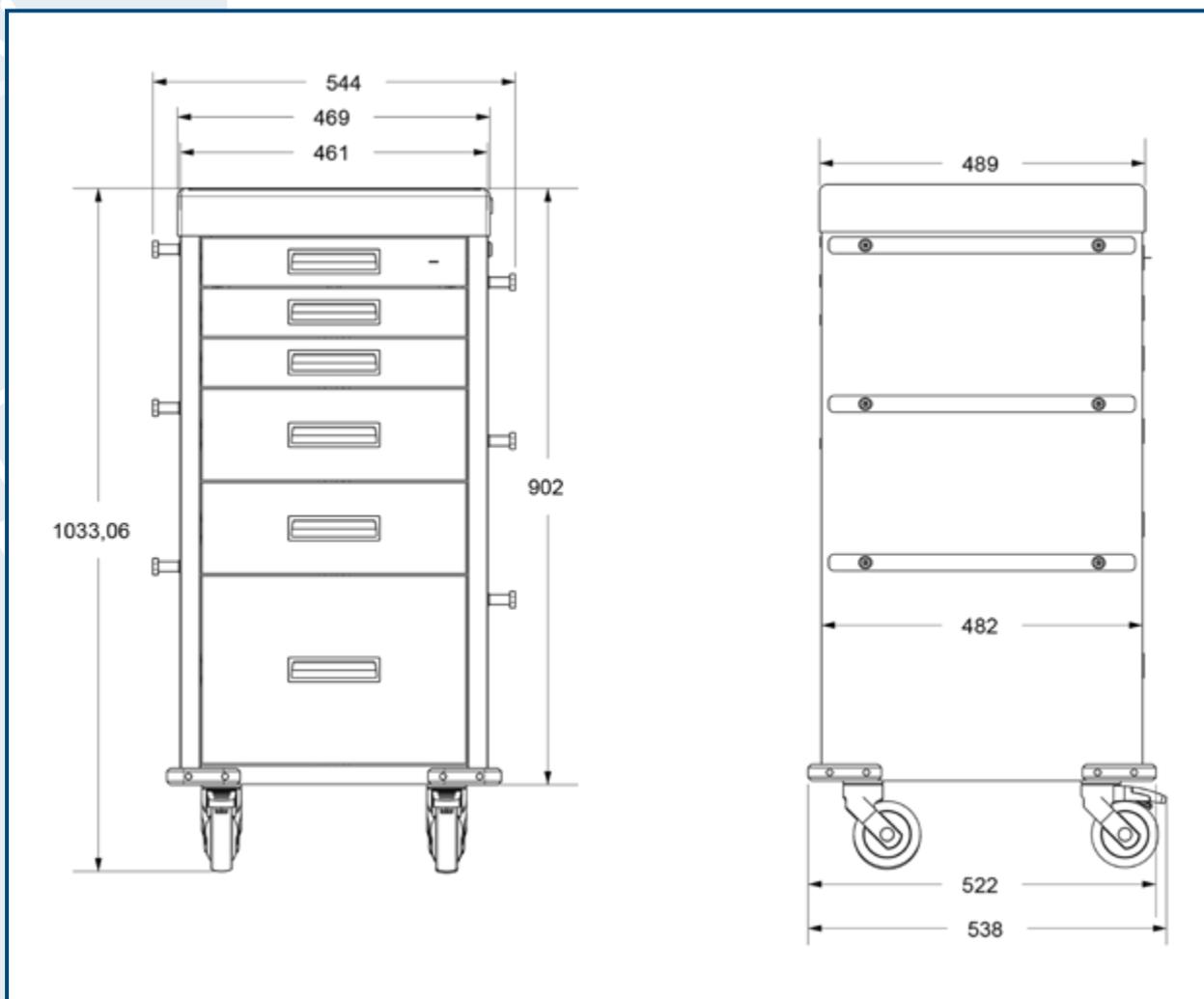

Technical sheet

Codice GMDN	46978
ID Code	GB0321.SP-1 GB0321.SP-2 GB0321.SP-3
Intended use	Multifunctional cart
Manufacturer	GARDHEN BILANCE SRL
Configuration 3+2+1 (standard)	3 drawers 42x37.8x7.8h cm 2 drawers 42x37.8x15.8h cm 1 drawer 42x37.8x23.7h cm
Configuration 1+3+1	1 drawer 42x37.8x7.8h cm 3 drawers 42x37.8x15.8h cm 1 drawer 42x37.8x23.7h cm
Configuration 2+1+2	2 drawers 42x37.8x7.8h cm 1 drawer 42x37.8x15.8h cm 2 drawers 42x37.8x23.7h cm

For certifications please refer to "Certifications and Quality system" section

Pictures contained in this brochure is for illustrative purpose only. The manufacturer reserves the right to modify the projects without notice. For accessories identification please refer to the accessory catalogue or commercial offer. When data are different please refer to commercial offer.

Dep. EN_189 rev. 28/06/2022

